



Web-based Software (ICDAS 3.0)

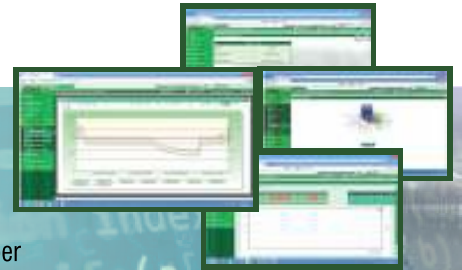

Access Anywhere Anytime

NEWTRONIC's ICDAS software changes the landscape of Stability Operations & Chamber Management. While it very effectively complies to the requirements of 21 CFR Part 11, it also meets the increasing demands for improved operational efficiencies. This server-based software gives the flexibility to operate the chambers remotely through any web browser and gives additional capabilities of going Paperless and system-based alert notifications of pending tasks.

Software Architecture

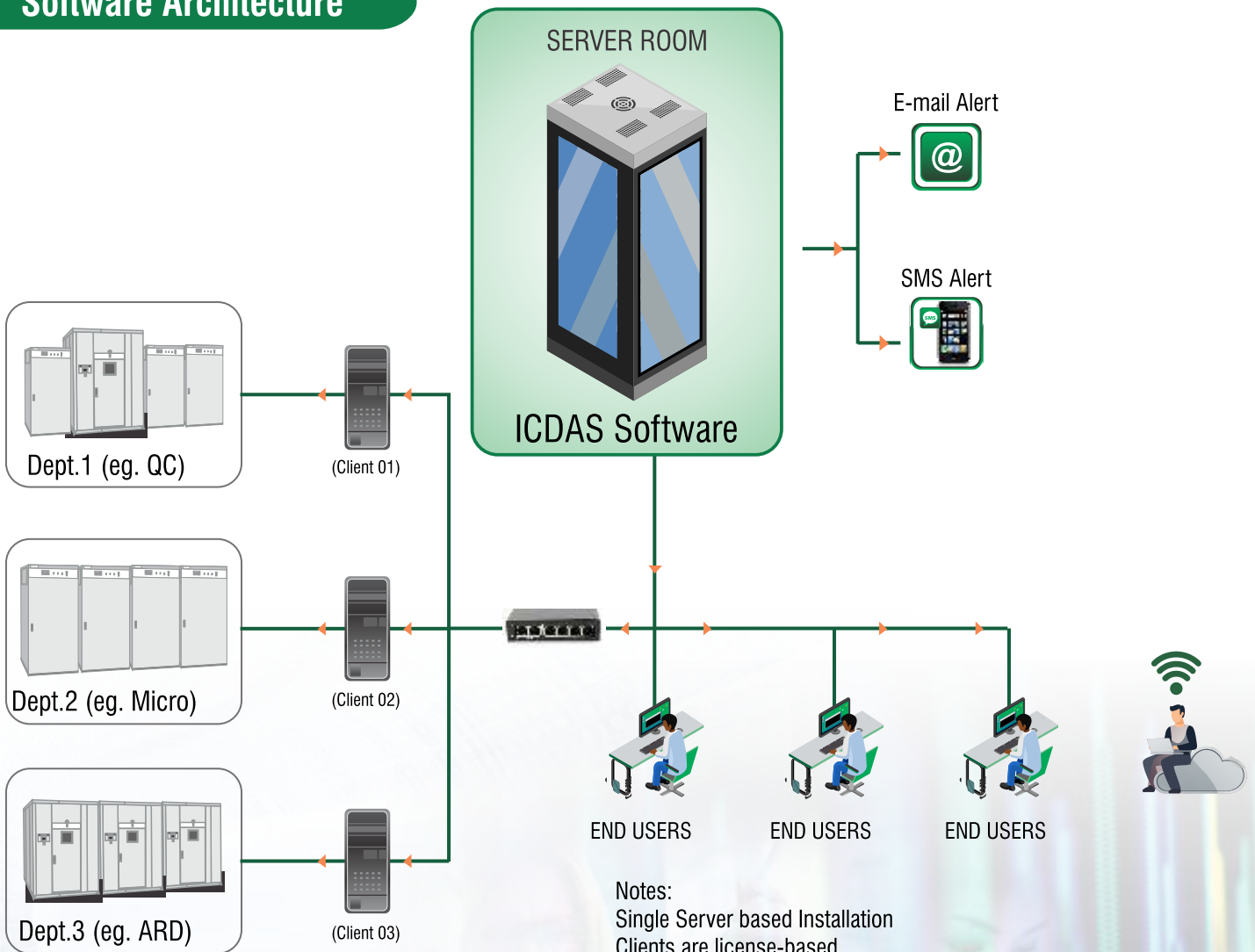

Notes:
 Single Server based Installation
 Clients are license-based
 1 client per dept. for upto 15 equipments
 ICDAS lite version available (in lieu of server based)



Ever-ready for surprise audits



Go Paperless



Restricted Access Rights & Privileges



Seamless Configuration




Safe & Secure Data



Date & Timestamp Audit Trails

Fast | Simple | Economical | Reliable

ALCOA+

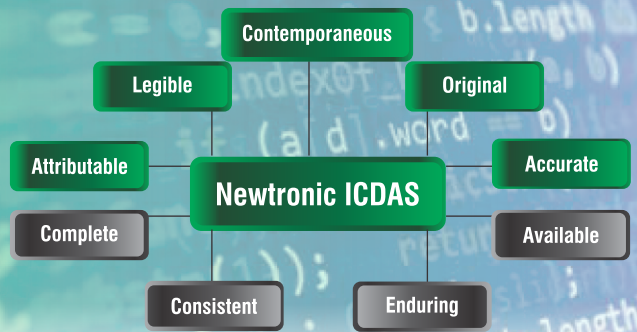

Instantly connect and control your chambers from anywhere without installing software.

Key Features



Installation

- Cloud-based Software accessible through Browser
- **Installation that is one-time for Lifetime!** Eliminate the need for multiple software installation & qualifications!



Data Management

- **GO PAPERLESS!!** e-Review, e-Approval, e-Signature, e-Comments of:
Data Logs | Alarms | Audit Trails
- Multi-level data authorization (through e-review & e-approval)
- Comprehensive Audit Trails of Alarms, User activity, Equipment settings, email & SMS generation
- **Stay Audit-Ready every single day through Supervisory Email facility!**
- Run-time calculation of MKT
- Real-time display of all connected equipment in single screen



Configuration

- Configure multiple user-groups (Operator, Supervisor, Admin, etc)
- Department-wise segregation of users, equipment & data
- Plant-wide harmonization of reporting format (header, footer, etc.)
- Configure the type & no. of alarms to be sent via email/SMS & configure delay time too



Safe & Secure

- Automatic & Manual daily data backup facility. Automatic data backup can be done both, remotely & locally
- Additional features of Data Archival & Data Restore



Other Features

- Viewing I/O Status & Running Hours of all major components. Component running hours to determine the shelf-life of major components
- Centralized user-defined Alarms, Alerts & Notifications through SMS & emails
- System Diagnosis for troubleshooting via Team Viewer
- Software to keep you USFDA 21 CFR Part 11 Ready



LDAP

- **LDAP (Light-weight Directory Access Protocol) Integration: Access ICDAS through your corporate email id/password**